

Probing Individual Proteins in Unsupported Membranes

Stephen W. Carmichael¹

Mayo Clinic

carmichael.stephen@mayo.edu

Proteins in biologic membranes perform a large variety of essential functions. The fact that about one third of all genes code for membrane proteins, and that the majority of drugs target these proteins, attest to that fact. However, until now, proteins have been studied under artificial conditions, such as after being crystallized, frozen, or adsorbed to a substrate. Rui Pedro Gonçalves, Guillaume Agnus, Pierre Sens, Christine Houssin, Bernard Bartenlian, and Simon Scheuring have devised a novel setup with the atomic force microscope (AFM) to allow proteins to be probed while they are in unsupported membranes.² Their method is similar in principle to methods where a small area of a membrane is sampled, such as when a piece of membrane is examined by patch clamping. The difference is that instead of attaching a membrane to the end of a pipette, it is spread across a piece of perforated Si(001). This provided "viewing holes" in the range of 90 nm to 250 nm. With the AFM, the membranes could be imaged and manipulated while in an aqueous buffered medium at ambient pressure and temperature.

Gonçalves *et al.* chose the surface membrane of *Corynebacterium glutamicum* as their model for imaging studies because this membrane has been so well characterized in many different ways. Patches of membranes could be clearly seen to be spanning holes on the nano-patterned Si(001) support, yet proteins within patches were considerably smaller than the holes. Using a specific AFM mode, they visualized the subunit architecture of the proteins. Such structural details as a pore about 15 Å in diameter,

corresponding to a channel pore, could be clearly seen. In addition to imaging, force measurements could also be made on the unsupported membranes. Specifically, they could analyze forces as the AFM tip made contact with the membrane and when it pierced the membrane. Measuring forces between these two events provided interesting information on the lateral interaction force between proteins within the membrane that suggested fairly weak protein-protein interactions reinforcing the importance of cooperative interactions.

Using *Halobacterium salinarium* purple membranes, known to contain proton pumps, Gonçalves *et al.* isolated wells of about 10 attoliters beneath suspended membranes. These wells were primed with the pH-sensitive fluorescent probe pyranine. Using quantitative fluorescence microscopy, they were able to demonstrate a significant increase in pH under appropriate circumstances. This was clear evidence that the pumping of protons was being directly observed.

This novel technique is an extension of several methods for imaging and manipulating individual membrane proteins as reviewed by Gonçalves and Scheuring.³ The key to their new method is that a membrane is suspended over a very small well. As techniques to manipulate fluids on a nano-scale are devised to mimic events occurring within cells, new methods can be developed to measure these events. Also, as AFMs become faster and more sensitive at acquiring images, these observations can be made at low physiologic levels and with better spatial and temporal resolution. We will be able to, in effect, see and measure what membrane proteins realistically do within cells with more precision than has ever been possible! ■

1. The author gratefully acknowledges Dr. Simon Scheuring for reviewing this article.

2. Gonçalves, R.P., G. Agnus, P. Sens, C. Houssin, B. Bartenlian, and S. Scheuring, Two-chamber AFM: probing membrane proteins separating two aqueous compartments, *Nature Methods* 3:1007-1012, 2006.

3. Gonçalves, R.P., and S. Scheuring, Manipulating and imaging individual membrane proteins by AFM, *Surface and Interface Analysis* 38:1413-1418, 2006.

INDEX OF ARTICLES

Probing Individual Proteins in Unsupported Membranes3

Stephen W. Carmichael, Mayo Clinic, Rochester, MN

Extraterrestrial Analyses of the Most Challenging Extraterrestrial Samples Ever Returned6

Frans J.M. Rietmeijer, U. of New Mexico, Albuquerque, NM

Microscopy Today 2007 Salary Survey Results12

Ron Anderson and Barbara Foster,* Microscopy Today and *Microscopy/Microscopy & Education, Inc.

Fibrous Sepiolite Use As An Asbestos Substitute: Analytical Basics18

Louis Solebello, Int'l Asbestos Testing Lab. Inc., Mt. Laurel, NJ

Low Vacuum SEMs: Latest Generation Technologies and Applications20

William Neijssen, Ben Lich, & Pete Carleson* FEI Company, Eindhoven, The Netherlands and *Hillsboro, OR, USA

Luminance Contrast – a New Visible Light Technique for Examining Transparent Specimens26

Jörg Piper, Clinic Meduna, Bad Bertich, Germany

Temporal Evolution of Incipient Damage on Metallic Surfaces Analyzed by Unidirectional Laser Oblique Illumination (ULOI)36

E.A. Favret^{1,2}, N.O. Fuentes^{2,3}, & L. Ferrero^{2,1}; Inst. de Tecnología Agropecuaria (INTA), ²Univ. Nacional de Gral San Martín, ³Comisión Nacional de Energía Atómica, Buenos Aires, Argentina

Multimode Trans-Illuminator for the Stereomicroscope ... 40

Ted Clarke, Retired Materials Engineer

New Approaches to Managing, Marketing, and Money for Maintaining a Core Facility

Developing a Financial Plan for the Long-term "Care and Feeding" of Major Equipment 44

Alberto Rodriguez and Charlene Sullivan, Krannert School of Management, Purdue University

The College of Microscopy – Meeting Rapidly Growing Microscopy Demands51

Donald A. Brooks, The McCrone Group, Westmont, Illinois

Getting Epoxy Semi-Thin Sections to Stick to Glass Slides52

Gilbert (Gib) Ahlstrand, Univ. of Minnesota, St. Paul, MN

SEM Stub Holders for Sputter Coating at 90° Tilt52

Amanda Best, Central Michigan University, Mt. Pleasant, MI

Digital Cameras and the TEM52

Warren Straszheim, Iowa State University, Ames, Iowa

Silicon Cross-Section Sample Preparation (Cleaving)53

Richard Beanland, Materials Analysis, Bookham, U.K.

Industry News 54

NetNotes57

Advertiser's Index78

ABOUT THE COVER

This composite image was taken by the navigation camera during the close approach phase of Stardust's Jan 2, 2004 flyby of comet Wild 2. Several large depressed regions can be seen. Comet Wild 2 is about five kilometers (3.1 miles) in diameter. To create this image, a short exposure image showing tremendous surface detail was overlain on a long exposure image taken just 10 seconds later showing jets. Together, the images show an intensely active surface, jetting dust and gas streams into space and leaving a trail millions of kilometers long. See the article by Rietmeijer starting on page 6. Image credit: NASA/JPL-Caltech.

Interested readers may want to Google "wild 2 comet pictures" to find many more images of Wild 2, including stereo pairs.

COMING EVENTS

2007

- ✓ **Microscopy and Microanalysis 2007**
August 5-9, 2007, Fort Lauderdale, FL
mm2007.microscopy.org
- ✓ **Cryo Technique and Immunogold labelling Workshop**
August 12-16, 2007, Emory University, Atlanta, GA
hyi@emory.edu
- ✓ **EMAG 2007 Electron Microscopy and Analysis Conference 2007**
September 3-7, 2007, Glasgow, Scotland
www.iop.org/Conferences... (search "forthcoming" then date)
- ✓ **EUROMAT 2007**
September 10-13, 2007, Nürnberg, Germany
www.euromat2007.fems.org/
- ✓ **Course: Live-cell imaging, FRET, FRAP and Spectral Unmixing**
September 18-21, 2007, Univ. of York, UK
www.york.ac.uk/depts/biol/tf/imaging_course.htm
- ✓ **FEMMS 2007 - Frontiers of E.M. in Material Science**
September 23-27, 2007, Sonoma, CA
femms2007.llnl.gov/index.html
- ✓ **Optical Microscopy and Imaging in the Biomedical Sciences**
October 9-18, 2007, MBL, Woods Hole, MA
www.mbl.edu/education/courses/special_topics/om.html
- ✓ **AVS 54th Symposium and Exhibition, (In-situ e.m. symp.)**
October 14-19, 2007, Seattle, WA
www.avs.org
- ✓ **Neuroscience 2007**
November 3-7, 2007, San Diego, CA
www.sfn.org/index.cfm
- ✓ **The American Society for Cell Biology**
December 1-5, 2007, Washington, DC
www.ascb.org

2008

- ✓ **ACCM-20 & IUMAS-IV**
February 10-15, 2007, Perth, Australia
microscopy.org.au/ACCM20/
- ✓ **PITTCON 2008**
March 3-6, 2008, New Orleans, LA
www.pittcon.org
- ✓ **American Soc. for Biochemistry and Molecular Biology**
April 3-9, 2008, San Diego, CA
www.asbmb.org
- ✓ **Scanning 2008**
April 14-16, 2008, Washington, DC
www.fams.org
- ✓ **Lehigh Microscopy School**
June 3-15, 2008, Bethlehem, PA
www.lehigh.edu/microscopy
- ✓ **Microscopy and Microanalysis 2008**
August 3-7, 2008, Albuquerque, NM
www.msa.microscopy.com
- ✓ **14th Electron Microscopy Congress, EMC 2008**
September 1-5, 2008, Aachen, Germany
www.euremicsoc.org/emc2008.html

2009

- ✓ **Microscopy and Microanalysis 2009**
August 3-6, 2009, Baltimore, MD
www.msa.microscopy.com

Please check the "Calendar of Meetings and Courses" in the MSA journal "Microscopy and Microanalysis" for more details and a much larger listing of meetings and courses.

MICROSCOPY TODAY

The objective of this publication is to provide material of interest and value to working microscopists!

The publication is owned by the Microscopy Society of America (MSA) and is produced six times each year in odd months, alternating with MSA's peer-reviewed, scientific journal *Microscopy and Microanalysis*. We greatly appreciate article and material contributions from our readers—"users" as well as manufacturers/suppliers. The only criterion is that the subject matter be of interest to a reasonable number of working microscopists. *Microscopy Today* has authors from many disparate fields in both biological and materials sciences, each field with its own standards. Therefore *MT* does not have a rigid set of style instructions and encourages authors to use their own style, asking only that the writing be clear, informative, and accurate. Length: typical article length is 1,500 to 2,000 words plus images, longer articles will be considered. Short notes are encouraged for our Microscopy 101 section. See our "Instructions to Authors" document on our website.

MICROSCOPY TODAY

ISSN 1551-9295

Ron Anderson, Editor

randerson20@tampabay.rr.com

Phil Oshel, Technical Editor

oshel1pe@cmich.edu

Thomas E. Phillips, Contributing Editor

PhillipsT@missouri.edu

Dale Anderson, Art Director

microscopytoday@tampabay.rr.com

Renée Stratmoen, Advertising Director

oshel1pe@cmich.edu

Regular Mail to:

Microscopy Today, P.O. Box 247, Largo, FL 33779

Courier Mail to:

1001 Starkey Road, Lot #374, Largo, FL 33771

Telephones:

1-(727)507-7101 • Fax: (727)507-7102 • Cell: (727) 631-1022

e-Mail:

microscopytoday@tampabay.rr.com

www Page:

<http://www.microscopy-today.com>



Colophon: *Microscopy Today* is created using components of Adobe Creative Suite CS2*

Total Circulation: 14,923

Disclaimer: By submitting a manuscript to *Microscopy Today*, the author warrants that the article is original (or that the author has the right to use any material copyrighted by others). The use of trade names, trademarks, etc., does not imply that these names lack protection by relevant laws and regulations. *Microscopy Today*, the Microscopy Society of America, and any other societies stated, cannot be held responsible for opinions, errors, or for any consequences arising from the use of information contained in *Microscopy Today*. The appearance of advertising in *Microscopy Today* does not constitute an endorsement or approval by the Microscopy Society of America of the quality or value of the products advertised or any of the claims, data, conclusions, recommendations, procedures, results or any information found in the advertisements. While the contents of this magazine are believed to be accurate at press time, neither the Microscopy Society of America, the editors, nor the authors can accept legal responsibility for errors or omissions.

© Copyright, 2007, The Microscopy Society of America. All rights reserved.

Bruker AXS Microanalysis

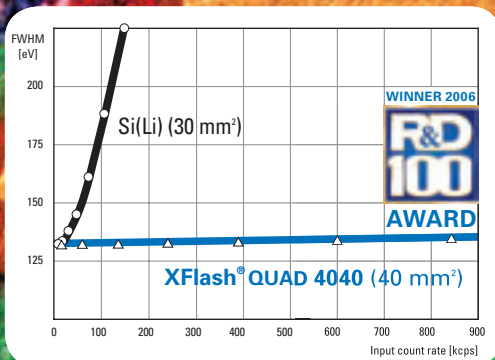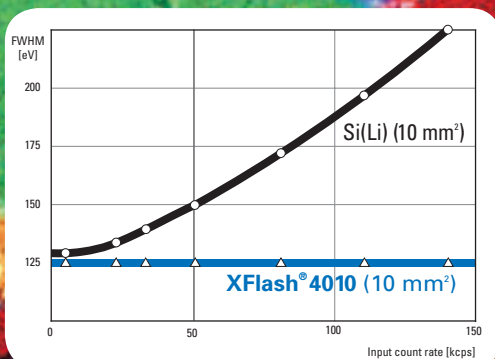

New
125 eV

Which EDS system is best for your microanalysis demands?

- 4th generation XFlash® silicon drift detectors (10, 30 and 40 mm²)
- Unmatched acquisition speed
- Excellent results at both low and high count rates
- Precise light element analysis (FWHM C-K α = 48 eV)
- New ESPRIT HyperMapping software (high speed PTS)
- LN₂ and vibration free

**QUANTAX – fast, reliable
and convenient microanalysis**

Visit us at the M&M 2007
booth # 320
to find out more about this
exciting new technology!

www.bruker-axs-microanalysis.com

think forward

X-RAY MICROANALYSIS