

Eucentric Stage for Recording Stereo Pair Photomicrographs

Ted Clarke


Metallurgical Failure Analysis Consultant

tclarke@ligtel.com

The stereo microscope is the most commonly used microscope. This applies to both the materials science and biological science fields. My past career as a materials engineer depended critically upon the three dimensional imaging of the stereo microscope for studying metal parts and fracture surfaces. Unfortunately our clients had no readily available way to view stereo pairs, so stereo

modernmicroscopy.com/main.asp?article=60. A simple tilting stage has traditionally been used by scientific photographers. This was not acceptable to me because the region of interest typically moves vertically and laterally when the stage is tilted. It is time-consuming to correct for these unwanted motions, especially when stereo pair recording is not frequently done. My solution was to build the eucentric stage shown in Figure 1 (left). The centering gage with a scribed line, shown on the lower right, is fitted to the spindle end of a water pump bearing. The camera is focused on the line and the stage is translated with the X-Y slides of the macro stand and rotated until the image of the scribed line bisects the field of the camera in the Y direction and is in sharp focus. The specimen


is then placed on the stage and the region of interest is brought into focus by using the vertical and horizontal slides and the rotational axis of the eucentric stage. Figure 2 (right) shows the back of the stage with the protractor and adjustable pointer. Figure 3 shows the first stereo pair recorded with the new stage with a tilt of 12 degrees between views. This stereo pair has been printed with a 60 mm center to center distance so it can be viewed with the more common viewers, such as the one included with the eighth edition ASM Metals Handbook "Fractography and Atlas of Fractographs." I

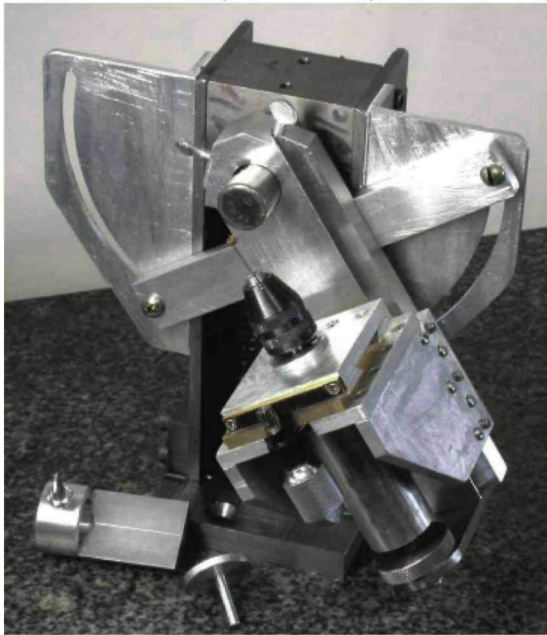

pair recording was not done. Likewise, our electron microscopes were not equipped with costly eucentric stages, which make stereo pair recording easily done. Jan Hinsch has helped me with the historical background on stereo photomicrography. Carl Zeiss offered a Druener camera attachment for their Greenough-style stereo microscope in 1900. More recently, Carl Zeiss, Wild and Nikon have offered sliding monocular tubes for recording either ray path with microscopes using a common main objective. The most recent technology, marketed by Buntun Instrument Co, Inc (www.buntgrp.com) allows video stereo imaging using a PC with live digital images from specially modified stereo microscopes equipped with dual digital cameras. This system serves as a vital training tool for surgeons. It appears that the entertainment industry will have to bring 3D imaging to the mass market, especially with PC based systems, before the technical community can make widespread use of communicating with stereo pair images. I suspect this will be a rapid transition.

I prefer viewing the separated prints with a large format Geoscope viewer purchased from Berezin Stereo Photography Products at: <http://www.berezin.com/3d/viewers1.htm> I am willing to work with anyone interested in making the eucentric stage commercially available. ■


Images from my recently completed scanning light photomicrography system convinced me of the need for stereo pair recording capability for this system shown in: <http://www.>

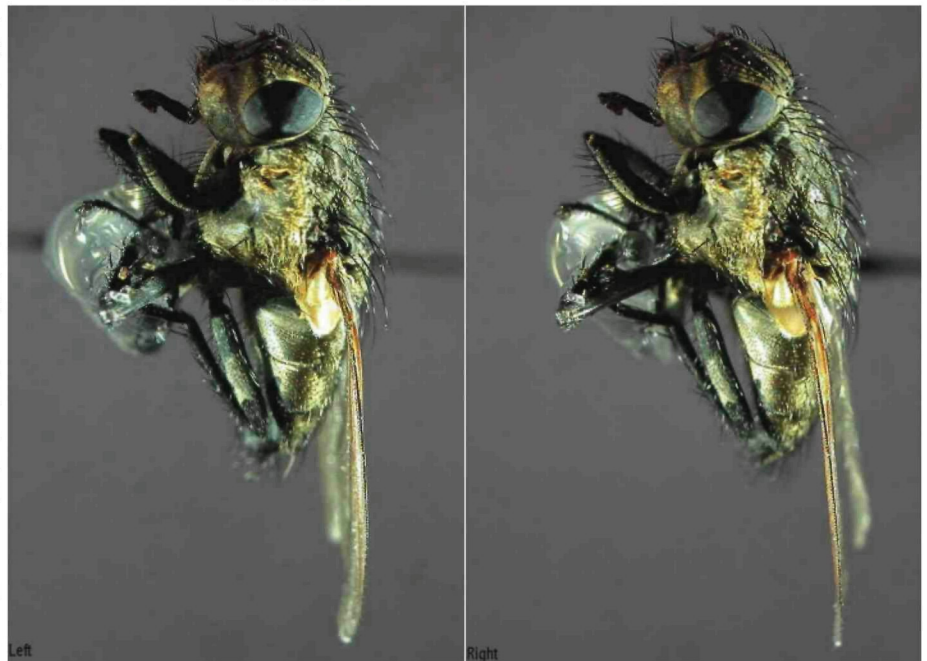

NEW


Omniprobe Short-Cut™ brings **HIGH THROUGHPUT** to in-situ TEM sample preparation

In-Situ FIB Lift-Out

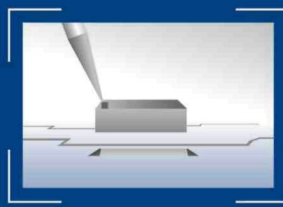

Streamlined in-situ FIB lift-out with Omniprobe's patented Total Release™ technology


Mechanical Conversion

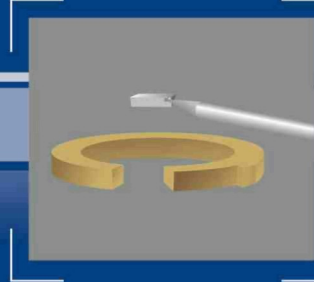


Direct mechanical conversion from lift-out to TEM sample

TEM sample produced using Omniprobe Short-Cut™ Tool

Omniprobe Short-Cut™ process for rapid TEM sample prep

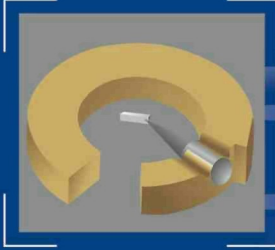

TEM Imaging

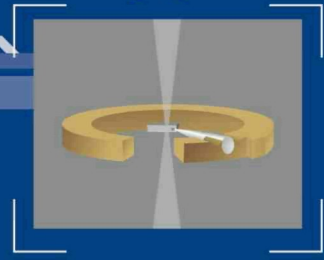

TEM Sample Prep And Nano-Mechanical Testing

The AutoProbe™ 200 is a computer controlled in-situ multi-purpose nanomanipulator system capable of:



- Lift-out for TEM Sample Prep
 - Electrical Measurements
 - Nano-Mechanical Testing
- Ideal For:
- IC Process Control and Failure Analyses
 - Topographic Surfaces
 - Brittle or Soft Materials
 - Phase Boundaries
 - MEMS Components

CE and UL certified


Patented Technology

AutoProbe™ 200

NEW

Computer Controlled Multi-Gas Injection System

The OmniGIS™ is Omniprobe's new multiple gas injection system for the FIB or SEM. The OmniGIS™ provides programmable control of three gas sources and an additional purge/carrier gas source.

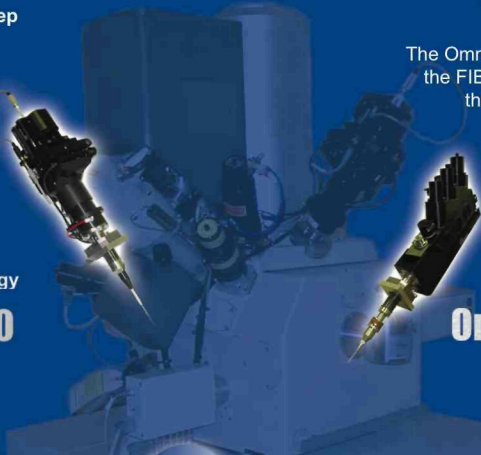

- Programmable Process Flow
- 3 Independent Gas Sources
- User-Replaceable Crucibles
- Feedback Control
- Patent Pending



Multiple Gas Sources (Oxide Isolated Pt Via)

OmniGIS™

Omniprobe, Inc.
10410 Miller Road
Dallas, TX 75238
Phone: 214-572-6800
Fax: 214-572-6801

omniprobe
Extending FIB Capabilities
www.omniprobe.com

European Support:
SEM-FIB Solutions
+44 1993 771203
tony.kendrick@sem-fib.com
Oxford, United Kingdom